M321318 V1

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DJO, LLC
1430 Decision Street
Vista
CA 92081
USA

©2011 DJO, LLC

Compex® WIRELESS

Quick Start Guide

1 DESCRIPTION OF THE DEVICE	2
2 FIRST USE: DEMO MODE	3
3 SYNCHRONISE YOUR REMOTE CONTROL	4
4 START WITH YOUR DEVICE	5
4.1 Select schedule, objective, setting or programme menu	6
4.2 Connect the modules to the electrodes	9
4.3 Preparing for your session	9
4.4 Start your session	10
5 RECHARGING	12

1 DESCRIPTION OF THE DEVICE

REMOTE CONTROL

- A** On/Off button (press briefly to switch on but press and hold >2s to switch off)
- B** 4 multi-function buttons:
- Functions associated with the icons shown on the screen (e.g. info, main menu, electrode placement, etc.)
 - Selection of stimulation channel to increase or decrease stimulation energy level
- C** Multi-directional pad (up-down-left-right) to select items
- D** Validation button
- E** Port for USB cable or docking station connector


MODULE

- A** On/Off button (press briefly to switch on but press and hold longer than 2 seconds to switch off)
- LED blinking green: ready to use
 - LED blinking yellow: on stimulation
- B** Groove to wind up the cable

DOCKING STATION

- A** Connector for charging the remote control
- B** Notch for opening the top of the docking station
- C** Slot for positioning the modules to be recharged
- D** Socket for the AC adaptor

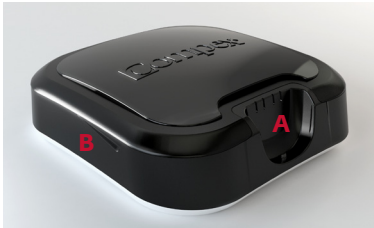2 FIRST USE: DEMO MODE

SELECT THE LANGUAGE

DEMO mode is available until you have paired your device with your account. DEMO mode only gives you access to certain programmes and not to all the possibilities of the Compex Wireless (empty calendar and no objectives).

Compex recommends that you visit www.compexwireless.com to create your own account which will allow you to download training plans and programmes, follow your progress on your calendar or create your own objectives.

Select your language

Validate your selection

DEMO MODE


demo Press DEMO button

3 SYNCHRONISE YOUR REMOTE CONTROL

CREATE YOUR PERSONAL ACCOUNT

In order to take advantage of all of the possibilities of the Compex Wireless (download training plans, create your own objectives, follow your progress, download new programmes, etc.), you should first create your personal account at www.compexwireless.com. Follow the instructions provided on the website. If you don't want to create an account, you can still use the product in DEMO mode (see item 2:First use: DEMO mode).

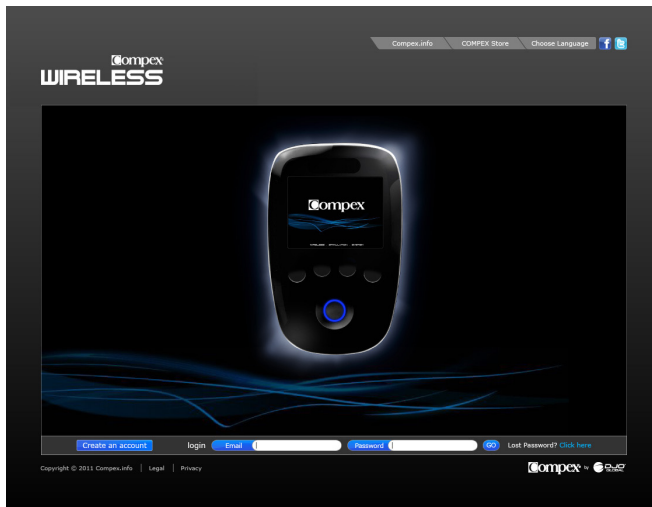PAIR YOUR REMOTE CONTROL WITH YOUR ACCOUNT

Once your personal account has been created, you can log into it and pair your remote control. To do so, please follow the instructions displayed on the website www.compexwireless.com

- A** Log into your account
- B** Download and install the software on your computer
- C** Click on the "Pair your remote" button and enter the serial number of your remote
- D** Connect the remote control to your computer using the USB cable provided and press the channel 4 button (the one with the icon ). A validation code will appear on the screen of the remote
- E** Type in the code in the designated field on the website
- F** Your remote control is now paired with your account

EN

4 GETTING STARTED WITH YOUR DEVICE

Once you have synchronised your remote control with your account and switched on your device, the first screen you see shows the tasks scheduled for that day. It is highly likely that you will see "day off" displayed on the screen, but if you had already scheduled Objectives, you would have seen some scheduled tasks displayed. From this screen you can begin your session (see Select an objective or a programme).

 **Back to main menu**

 **Press for Info menu**

Select an item

Validate your selection

 **Back to main menu**

4.1 SELECT SCHEDULE, OBJECTIVES, SETTINGS OR PROGRAMMES MENU

Select your item

Validate your selection

1. SCHEDULE

Here you can find your schedule with the session you did or you will have to do.

- A green rectangle below the day means that you have done what was planned
- An orange rectangle below the day means you have not done what was planned in the past
- A grey rectangle below the day means that you have planned an activity in the future

Select a day to view the details for it. From there you can begin your session (see Select an objective or a programme).

Validate your selection

Select a day

 Back to main menu

Validate your selection

Select an item

 Previous screen
 Press for Info menu

EN

2. OBJECTIVE

Here you will find the objectives you have downloaded from the website. To learn how to navigate in this menu, see Select an objective or a programme.

4. PROGRAMMES

The programmes menu contains the programmes you have downloaded from the website.

To learn how to navigate in this menu, see Select an objective or a programme. Please note: Only a limited number of programmes are installed in DEMO mode. Once you have paired your remote control with your account, you will be able to modify this list by adding or removing programmes.

3. SETTINGS

In the settings menu you can set several parameters such as backlight, volume, language, etc.

Validate your selection

Select your item

Back to main menu and apply changes requested

SELECT AN OBJECTIVE OR A PROGRAMME

The steps for selecting an objective or a programme are similar.

SELECT AN OBJECTIVE

Validate your selection

Select your item

Back to main menu

SELECT A PROGRAMME

Validate your selection

Select your item

Back to main menu

Press for Info menu

When an objective is selected you can see the details of all the sessions it comprises.

Back to previous screen


Select your item

Validate your selection

Press for Info menu

See the electrode placement and body position associated with your objective and selected task. For programmes, you will be able to view all the possible electrode placements. For objectives, you will see what placement to use. To scroll through the illustrations, use the up/down button on the multi-directional pad.

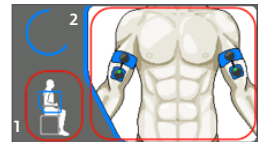You can now place the electrodes on your body and connect the modules to them.

Please note: For voluntary sessions (cycling in this example), you have to validate the session yourself by pressing the button under the icon

1. Body position
2. Electrode placement



EN

4.2 CONNECT THE MODULES TO THE ELECTRODES

Stick the electrodes on your skin. Electrodes are attached from the side. Slide the module along the electrode snap until you hear a click. Please note: Markings tell you which way to insert the On/Off button on the main pod and a little vertical line on the shell of the other pod.

To disconnect the module from the electrode, just do the reverse.

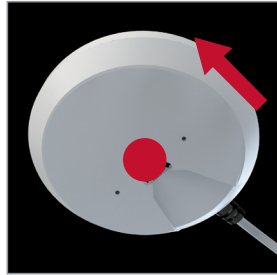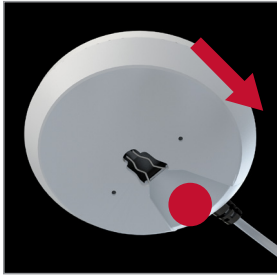4.3 PREPARING FOR YOUR SESSION

WARM-UP AND COOL DOWN

This screen is not available for all programmes. You can choose to do the warm-up phase before starting the session and the cool-down phase at the end of the stimulation session.

Please note: Compex recommends that you keep the warm-up phase, unless you have already done a voluntary warm-up, and the cool-down phase.

Select (or not) the warm-up or cool-down phase

Back to previous screen

Validate your selection

WIRELESS

4.4 START YOUR SESSION

SWITCH ON THE MODULES BEFORE STARTING THE STIMULATION SESSION

Before initiating the stimulation and mi-SCAN measurement, you have to turn on the modules one after another by pressing their On/Off buttons. Once the 1st module is detected by the remote control, you can switch on the 2nd module and so on up to the 4th module. If you only want to use two modules for your session, just press the START button after the second module is detected.

Please note: For a session which is optimally adapted to your muscles, Compex recommends that you always perform an mi-SCAN measurement before starting a stimulation session (see explanation below).

-  **start** Start your session without the mi-SCAN measurement
-  **mi** Start your session with the mi-SCAN measurement (see explanation below)
-  **◀** Back to previous screen

BEGINNING OF THE STIMULATION SESSION: MI-SCAN MEASUREMENT

This feature adapts the electrostimulation session to each individual user's physiology. Just before starting your stimulation session, mi-SCAN probes the chosen muscle group and automatically adjusts the stimulator parameters to the excitability of this area of the body. It is a truly personalized measurement. In order to ensure that your session is as effective and comfortable as possible, it is highly recommended that you perform the mi-SCAN at each session.

This function is implemented at the start of the programme by a short sequence during which measurements are taken. You should remain completely still and be perfectly relaxed throughout the measurement test. Once the test has been completed, the programme can start.

EN

DISPLAY EXPLANATIONS

Remaining time in minutes

Remaining number of contractions or indication of warm-up or cool-down phase

Baregraph of energy level

Energy level (goes up to 999)
Number on black background: active background
Number on grey background: inactive channel

PROGRAMMES INDUCING MUSCLE CONTRACTIONS
Top sector: contraction phase
Bottom sector: active rest phase

PROGRAMMES NOT INDUCING MUSCLE CONTRACTIONS
The top and bottom sectors show the time remaining in the specific phase

DURING THE STIMULATION SESSION

Stimulation always starts off at 000 (+++).

Pause the stimulation

Increase or decrease the stimulation energy levels of the selected channels

Select the channels on which you wish to act. When a channel is active, the LED of the button gives off a solid blue light.

Please note: by default, all channels are active at the start of a session.

Please note: pressing the button for a given channel for more than 1 second activates that channel and deactivates all the others.

266 Active channel (Black background)

266 Inactive channel (Grey background)

Please note: Active rest phase stimulation energies are automatically set at 50% of contraction energies. These can be modified during the rest phase. Once modified, they will be totally independent of the contraction energies.

PAUSE MODE

You can pause the device by pressing the centre button or the On/Off button on one of the modules during the stimulation.

Please note: The session will resume at 80 % of the energy levels that were being used prior to the interruption.

 **start** Resume the stimulation session

 Back to the screen displayed before the stimulation session

END OF SESSION

At the end of the session, a screen with a check mark comes up. Press any button to return to the main menu. To switch off the device, just press and hold the On/Off button on the remote control for more than 2 seconds. This will also switch off all the modules.

5 RECHARGING

BATTERY LEVEL

The battery level of each module is shown when you switch them on before starting your stimulation session. The battery level of the remote is always visible in the upper right corner. The little green circles indicate how many modules are switched on and can be detected by the remote control.

Remote control battery level
Number of modules which are on
and detected by the remote control
Module battery level

EN

CONNECT THE DOCKING STATION

Connect the AC adaptor supplied with your device to the docking station then plug it into a socket. Compex strongly advises you to fully charge the remote and module batteries before using them for the first time as this will improve their performance and life span.


CHARGING THE MODULES AND THE REMOTE CONTROL

At the end of your stimulation session, it is highly recommended that you store the remote control and the modules in the docking station in order to recharge them.

THIS SCREEN COMES UP AS SOON AS THE REMOTE CONTROL IS CONNECTED TO THE DOCKING STATION.

Then, place the modules in the spaces provided.

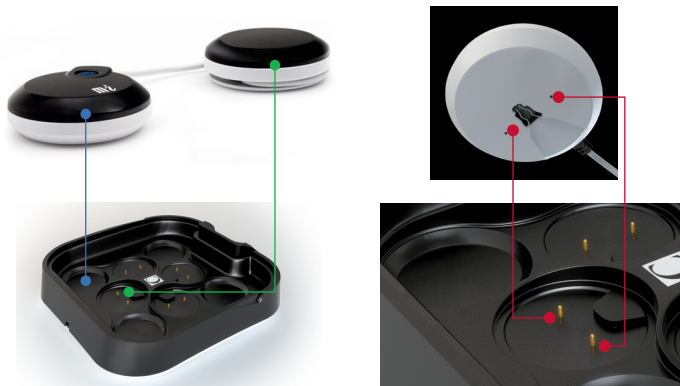To do this, place the pod without the On/Off button (the **green** one in the illustration) in the slot indicated in **green** and the other pod in the **blue** slot. Do the same for the others modules.

The pod without the On/Off button has to be inserted over the small connectors. A magnet and the little vertical mark on the shell of the pod should help you to properly position the pod in its slot. When placed correctly, you should hear a click.

When you place a module on the docking station, a red circle will appear to indicate that charging is in progress. The module is still charging as long as the red circle is visible. As soon as this disappears, the module is completely charged and ready for use.

Please note: If you place a fully charged module on the docking station, no red circle will appear.

As soon as the remote control and the modules are fully charged, the current time will be displayed.

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